




ENGINEERED ELECTRONICS, INC.

FEATURES

- * 80 POINT CAPACITY
- * DESK TOP, SURFACE, OR FLUSH MOUNTING
- * SIGNALS 4 LEVELS OF CALL
- * EMERGENCY CALL IDENTIFICATION
- * POWER ON INDICATION
- * SOFTWARE FOR ROOM LABELING
- * SPILL PROOF LAMP ASSEMBLY
- * LED LAMP TEST

NURSE CALL ANNUNCIATOR 900-510

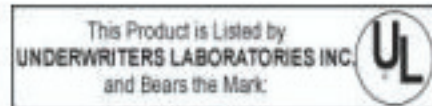

DS2000 V

GENERAL DESCRIPTION

The EEI DS2000 V Annunciator is designed to receive calls from Patient or Emergency Stations and indicate the level of call with a pulsing tone and either pulsing or steady light depending on type of call for individual room identification.

ENGINEERS' SPECIFICATIONS

- (1) All equipment described shall be manufactured and assembled entirely in the United States of America.
- (2) The Annunciator shall have capacity to accommodate 80 individual room positions. (30 positions - 900-513, 50 positions - 900-515)
- (3) Each room button shall be back lit with a red LED.
- (4) Lamp tests of all LEDs are accomplished from front of panel.
- (5) All room positions can be labeled at factory or through furnished software.
- (6) To prevent damage from spillage, the Annunciator shall be constructed so that drinks can not be set on the unit. If liquid does get on front panel, it shall not affect operation or get to any electronic components.
- (7) There shall be an emergency indicator lamp which will illuminate red whenever any emergency device is activated.
- (8) There shall be a green power on indicator.
- (9) The cabinet shall be constructed of heavy gage painted steel with a lexan overlay across the entire front. The cabinet shall be designed for desktop, surface, or flush mounting.
- (10) Plugging in supplied 25 foot interconnect cables shall be the only connection necessary to connect Annunciator to Central Control Panel.
- (11) The DS2000 V Annunciator shall conform to Underwriters' Laboratory, Inc. Standard (U.L. 1069).



Standard Equipment

One - 500-057 25 ft. Interconnect Cables (900-513)
Two - 500-057 25 ft. Interconnect Cables (900-510,900-515)

MODELS

| | |
|---------|----------------------|
| 900-513 | 30 point Annunciator |
| 900-515 | 50 point Annunciator |
| 900-510 | 80 point Annunciator |

Dimensions

| | |
|-------------------|---------------------------|
| Surface mounting | 11 3/4"H x 14"W x 3"D |
| Flush mounting | 11 3/4"H x 14"W |
| Desk top mounting | 11 3/4"H x 14"W x 8 1/4"D |

NOTES

- (1) One or two twenty-five foot 500-057 interconnect cables are supplied with each Annunciator as a standard feature. Longer cables are available upon request. The cables are supplied with connectors on each end.
- (2) A minimum of 1 1/4" conduit is required between Annunciator and Central Control Panel to accommodate the 500-057 interconnect cables.

FLUSH CUTOUT SIZE

